

The Prelude to your Symphony



Dia-PrepTM *Plus*




EDTA Root Canal Preparation Cream

DiaDent[®]

DiaDent Group International

#11-3871 North Fraser Way, Burnaby, BC Canada V5J 5G6 | T: 604.451.8851 |

TF: 1.877.DIA.DENT (1.877.342.3368) | F: 604.451.8865 | diadent@diadent.com | www.diadent.com

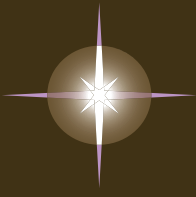

Just like a *prelude* to a symphony, Dia-Prep Plus prepares you for the success of your root canal treatment. DiaDent Dia-Prep Plus is an **EDTA** cream with Urea Peroxide, used as an irrigation agent for chemo-mechanical preparation of the root canal.

advantages & effectiveness

- EDTA removes smear layer, softens dentin, and facilitates the removal of calcific obstruction
- Urea peroxide reacts with sodium hypochlorite to produce oxygen bubbles that flush out debris within the canal
- The lubricant medium facilitates instrument movement within the canal to eliminate instrument binding or breakage
- Unclogged dentin tubules and apical area permit better penetration of root canal medicaments and fillings

order information

Item # 2002-1100

- 2 x 6 g Dia-Prep Plus syringes
- 10 disposable tips
- 1 protective cap


indications

(use in conjunction with sodium hypochlorite)

- Debridement in the root canal
- Dissolution of pulp remnants
- Removal of smear layer
- Facile movements of canal instruments
- Microbicidal treatment within the canal
- Clearance of dentin tubules and apical area

composition

- Ethylenediamine Tetraacetic Acid (EDTA) 15%
- Urea Peroxide 10%
- Polypropylene Glycol & Polyethylene Glycol

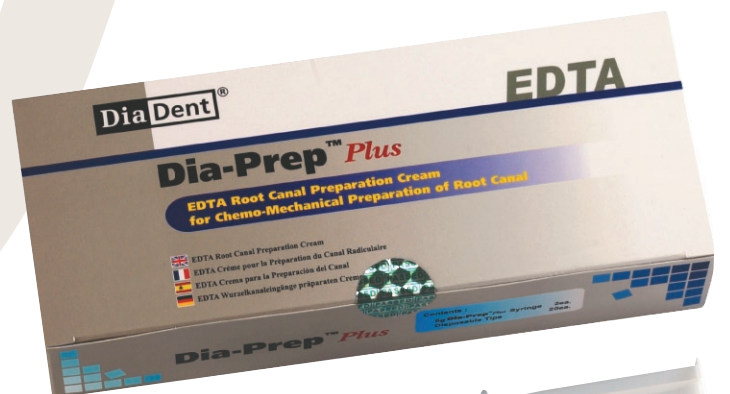

DiaDent®

www.diadent.com
1.877.342.3368